


TRIVA-aA /-aAY

TECHNICAL INFORMATION


TRIFOCAL
MONOBLOC
ACRYLIC
ASPHERIC
PRELOADED



TRIVA-aA

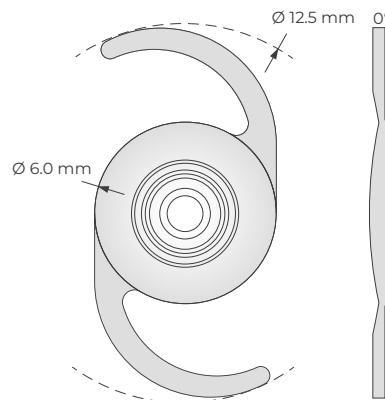

TRIVA-aAY
YELLOW

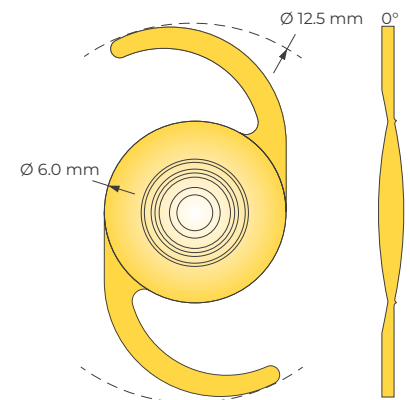

Type

- foldable, trifocal posterior chamber IOL
- one-piece with C-loop haptics
- available in preloaded Safeloader autoloading container or Compact Line

Material

- hydrophilic acrylic with UV-absorber
- water content 26 %
- blue-light protection optional
- Abbe number 56
- refractive index 1.46

Specifics

- central diffractive aspheric anterior surface with a refractive optic periphery
- intermediate addition + 1.75 D, near addition +3.5 D (at IOL plane)
- posterior surface with 360° lens epithelial cell barrier

Diopter Range

Sphere in 0.5 D steps
(optic shape)

10.0 to 30.0 D (biconvex) **SL** **CL**

Est. A-const

For a complete list of recommended constants,
please visit www.humanoptics.com



SL = Safeloader | CL = Compact Line

CE0044

 HumanOptics Holding AG

Spardorfer Str. 150 · 91054 Erlangen · Germany · Phone +49 (0) 9131 50 66 5-0 · Fax +49 (0) 9131 50 66 5-90
mail@humanoptics.com · www.humanoptics.com